

# Vision tests



- SZB Bailey-Lovie chart set
- Set of 2 panels with logarithmic progression of optotypes
- 5 optotypes per visual acuity level
- Suitable for use at any viewing distance
- Panel 1 for measuring distance visual acuity at normal contrast
- Panel 2 for measuring distance visual acuity at low contrast
- Dimensions: 56x65 cm (WxH)

**938147** SZB Bailey-Lovie chart set


**938147-1** SZB Bailey-Lovie chart set

each panel divided into upper and lower half



**938157-EN** Near vision reading test set

to determine the magnification requirement for near vision including Amsler grid test

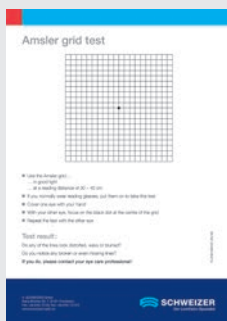

**938107-EN** Amsler grid test

## ■ Near vision reading test set

Up-to-date test in logarithmic layout according to DIN 58220 and ISO 8596. The test chart consists of unrelated sentences in different letter sizes. The notation below the smallest print the patient can still read fluently indicates the magnification required for reading newsprint. The Amsler grid test is included in the set.

## ■ Amsler grid test

The Amsler grid test for detecting and monitoring symptoms of macular degeneration is straightforward to administer. It comes printed on sturdy paper for ease of use.